

Computerized Tonometer

CT-800



 **TOPCON**[®]
CONNECTING VISIONS

CT-800

COMPUTERIZED TONOMETER

Smooth Operation

Topcon is proud to introduce the new CT-800 Computerized Tonometer which incorporates the very latest in design technology and ergonomics. The unit features a bright, new 8.5 inch color touch screen panel to control the main functions and an improved joystick operation due to a 22% reduction in weight and motorized up/down operation from the previous

model. This model provides a soft air-puff for measurement, which equates to less stress on the eye. Topcon is always striving to deliver the very best instrumentation which always incorporates the best design technology, functionality and accuracy, allowing our customers to have confidence in their investment for years to come.

FEATURES

- » Smooth operation
- » Compact & Modern Design
- » 8.5 inch LCD Touchscreen Panel
- » Soft Air Puff
- » Adjusted Tono Value*
- » Connectable LAN



COMPUTERIZED TONOMETER **CT-800**

SMOOTH OPERATION

Topcon's ability to engineer a weight reduction of approximately 22 % has contributed to a smoother operation of the unit when using the X-Y control lever. The up/down movement has also been motorized. This all contributes to a more fluid measuring process and a more accurate end result.



WIDE 8.5-INCH TOUCHSCREEN CONTROL PANEL

The new CT-800 is equipped with a wide 8.5 inch color touchscreen panel featuring clear, easy to read icons allowing the operator full control of the unit during the measurement process.



SOFT AIR PUFF

Software controls the optimal amount of air required for every measurement for each individual eye. A soft air puff means less stress on the eye, resulting in an improved patient experience.

ADJUSTED TONO VALUE*

The CT-800 is able to calculate Adjusted Tono Value from measured IOP to manually input the corneal thickness value measured by other devices such as a pachymeter.



* Adjusted Tono Value is reference value

IOL MODE

The IOL mode allows you to measure the eye with an IOL.

ONE TOUCH LOCK, EASY LOAD PRINTER AND NOZZLE COVER


The CT-800 is equipped with a built-in thermal printer that can be easily loaded with a roll of paper. The unit automatically cuts the paper at the end of printing. And the nozzle cover can be simply attached to the CT-800 body.



Specifications

Ocular Pressure Measurement	1 to 60 mmHg (1 mm Hg step)
Measuring Range	Average indication 1mmHg / 0.1mmHg steps changeable
Measuring Mode	1-30 mmHg / 1-60 mmHg
Monitor	8.5" WVGA color LCD monitor, touch panel
Data Transport Terminal	USB (Import), RS232C (Export), LAN (Export)
Dimensions	317 - 341mm(W) × 521 - 538mm(D) × 437 - 467mm(H)
Weight	14.0 kg
Power Supply	100-240V AC, 50-60Hz, 30-70VA

System chart

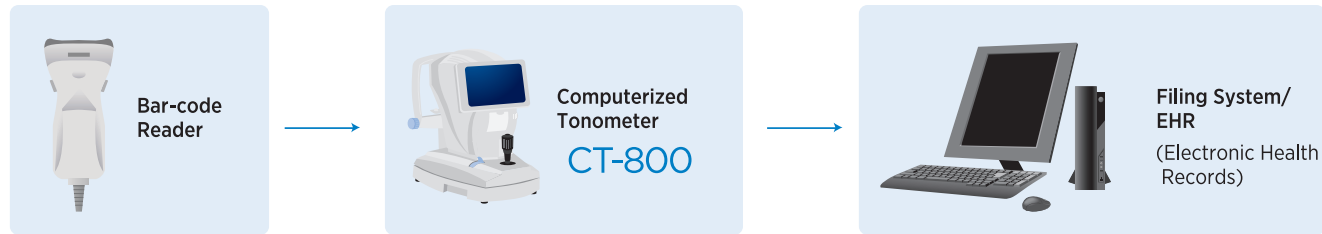


* Contact Topcon subsidiaries or dealers for system configurations.

Printed Data

```

    2014_02_14      AM 10:37
    TONO. DATA
    <R>             <L>
    mmHg  ADJ.  mmHg  ADJ.
    17    17    17    17
    16    16    16    16
    17    17    17    17
    AVG.
    16.7  16.7  16.7  16.7
    PACH. DATA
    <R>  mm      <L>  mm
    0.535      0.532
    TOPCON
  
```

- 1 Tono value on display unit "mmHg"
- 2 Tono average value on display unit "mmHg"
- 3 Adjusted Tono value
- 4 Pacho value with manual input

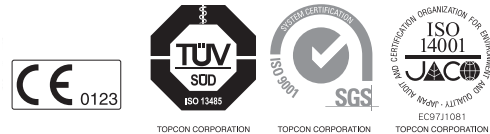

IMPORTANT Subject to change in design and/or specifications without advanced notice.
 In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation.

Topcon Europe Medical B.V.
 Essebaan 11; 2908 LJ Capelle a/d IJssel; P.O. Box 145;
 2900 AC Capelle a/d IJssel; The Netherlands
 Phone: +31-(0)10-4585077; Fax: +31-(0)10-4585045
 E-mail: medical@topcon.eu; www.topcon-medical.eu

Topcon Denmark
 Praestemarksvej 25; 4000 Roskilde, Danmark
 Phone: +45-46-327500; Fax: +45-46-327555
 E-mail: info@topcondanmark.dk
 www.topcondanmark.dk

Topcon Scandinavia A.B.
 Neongatan 2; P.O. Box 25; 43151 Mölndal, Sweden
 Phone: +46-(0)31-7109200; Fax: +46-(0)31-7109249
 E-mail: medical@topcon.se; www.topcon.se

Topcon España S.A.
 HEAD OFFICE; Frederic Mompou, 4;
 08960 Sant Just Desvern; Barcelona, Spain
 Phone: +34-93-4734057; Fax: +34-93-4733932
 E-mail: medica@topcon.es; www.topcon.es

Topcon Italy
 Viale dell' Industria 60;
 20037 Paderno Dugnano, (MI) Italy
 Phone: +39-02-9186671; Fax: +39-02-91081091
 E-mail: info@topcon.it; www.topcon.it

Topcon France
 BAT A1; 3 route de la révolte, 93206 Saint Denis Cedex
 Phone: +33-(0)1-49212323; Fax: +33-(0)1-49212324
 E-mail: topcon@topcon.fr; www.topcon.fr

Topcon Deutschland GmbH
 Hanns-Martin-Schleyer Strasse 41;
 D-47877 Willich, Germany
 Phone: (+49) 2154-885-0; Fax: (+49) 2154-885-177
 E-mail: med@topcon.de; www.topcon.de

Topcon Polska Sp. z o.o.
 ul. Warszawska 23; 42-470 Siewierz; Poland
 Phone: +48-(0)32-670-50-45; Fax: +48-(0)32-671-34-05
 www.topcon-polska.pl

Topcon (Great Britain) Ltd.
 Topcon House; Kennet Side; Bone Lane; Newbury
 Berkshire RG14 5PX; United Kingdom
 Phone: +44-(0)1635-551120; Fax: +44-(0)1635-551170
 E-mail: medical@topcon.co.uk; www.topcon.co.uk

Topcon Ireland
 Unit 276, Blanchardstown; Corporate Park 2
 Ballycoolin; Dublin 15, Ireland
 Phone: +353-18975900; Fax: +353-18293915
 E-mail: medical@topcon.ie; www.topcon.ie



TOPCON CORPORATION
 75-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8580, Japan.
 Phone: 3-3558-2523/2522. Fax: 3-3960-4214. www.topcon.co.jp

Item code: 5213591 / printed in Europe 0314