

GoVu™


PACHYMETER



Hand-held and easy to use

Instantly and accurately measures
central corneal thickness

Battery-operated with
charging cradle included

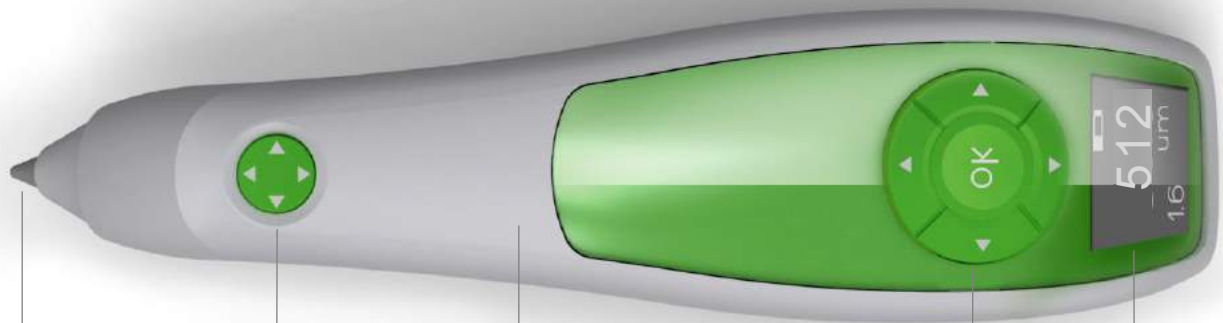

Integrated 20MHz transducer for excellent accuracy

Full scan controls available at fingertips

Ergonomic shape comfortably fits either hand

Easy navigation of scan results and features

Large screen for viewing results and battery status

Charge. Then Go Vu Central Corneal Thickness.

Specifications	Length	7.5"
	Weight	3.6 oz
	Display	1.1" Color LCD
	Battery	Rechargeable Lithium-Ion 3.7V / 750 mAh
	Charger	Input: 90-264 VAC (Universal) Output: 5 VDC / 1A
	Printer	Optional Bluetooth thermal printer
	Transducer	20 MHz frequency
	Measurement Range	300 - 1000 μm
	Clinical Accuracy	Within 5 μm
	Measurement Clock	160 MHz
Averaging	25 measurements per displayed result	
Precision	1 μm	