

## **IRIDEX®** LIO Plus

An Easy to Use Single-Mirror Laser Indirect Ophthalmoscope
That Is Ergonomic and Provides Excellent Peripheral Visualization,
Treatment Flexibility, Consistency, and Reliability



For many years, the IRIDEX® LIO Plus, a single-mirror laser indirect ophthalmoscope, has been considered **the gold standard of LIOs**. Retina specialists and teaching institutions around the world use it for its excellent visualization, ease of use, and reliability.

- The IRIDEX LIO Plus facilitates observation and photocoagulation of peripheral retinal pathologies including diabetic retinopathy and retinopathy of prematurity<sup>1</sup>
- Provides excellent peripheral visualization and consistent laser uptake as far as the ora serrata
- It is ergonomic, intuitive, and ideal to use in patients who are best examined and treated in a supine position
- Has a lightweight headset with size adjustments to optimize fit and comfort
- Permits independent positioning of the laser spot within the illuminated field, or simultaneous adjustment of both the laser spot and illuminated field within the user's visual field
- The long laser depth of focus tolerates a wide range of working distances
- The 810 nm Large Spot model can be used with the OcuLight SLx Laser to deliver transpupillary thermotherapy for the treatment of intraocular tumors.<sup>2,3</sup>
- The single-wavelength models are compatible with 532 nm, 577 nm, and 810 nm Iridex lasers
- o A dual-wavelength model can be used with both 532 nm and 810 nm Iridex lasers for added utility
- The halogen illumination source offers consistent illumination and great color rendering
- Its integrated laser filters provide protection for the physician during laser use
- Includes a travel case for safe storage and ease of transportation to satellite offices




"Treating comfortably, with minimal ergonomic burden, is an important way to help ensure


consistency of results across treatments. The LIO Plus, with its intuitive design and practical features,

plays a large role in my laser treatment protocol — It is proving to be effective for the most challenging

peripheral cases, and my colleagues and I have been pleased with the results."

Shawn Kavoussi, MD
Texas Retina Center, Houston, United States

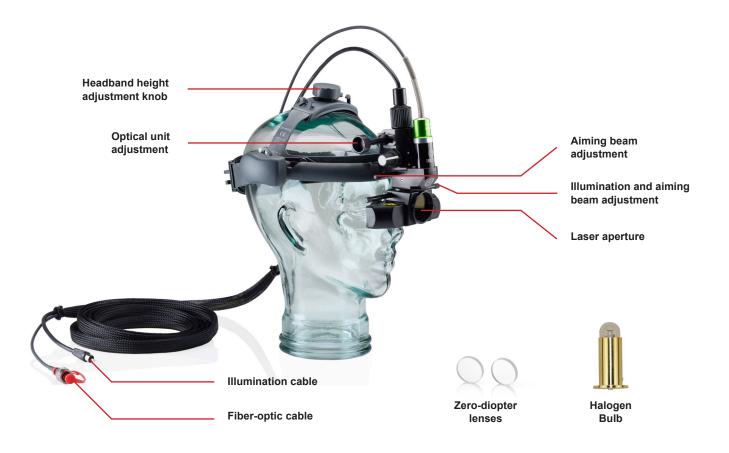

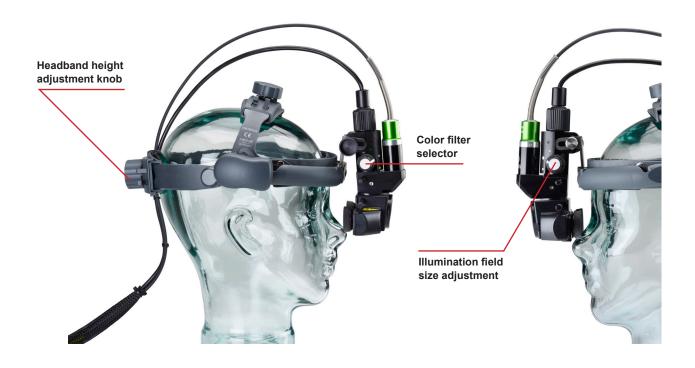

## Compatibility, Specifications, and Ordering Information

| Product Number            | Wavelength      | Spot Size* | Laser Compatibility                   | Connector |
|---------------------------|-----------------|------------|---------------------------------------|-----------|
| 30903 - H500              | 532 nm / 810 nm | 360 μm     | OcuLight GL, GLx, TX, SL, SLx, IQ 532 | Resistive |
| 13152 - H500              | 810 nm          | 360 µm     | OcuLight SL, SLx, DioVet              | Resistive |
| 13153 - H500 (Large Spot) | 810 nm          | 1400 μm    | OcuLight SL, SLx, DioVet              | Resistive |
| 65515 - H500              | 532 nm          | 360 μm     | IQ 532                                | RFID      |
| 65900 - H500              | 577 nm          | 360 µm     | IQ 577                                | RFID      |

| Dimensions |                   | Ophthalmoscope                          |                        |  |
|------------|-------------------|---|------------------------|--|
| Width      | 7.0 cm / 2.75 in  | Headband circumference adjustment range | 520 mm to 640 mm       |  |
| Length     | 13.5 cm / 5.31 in | Illumination field sizes**              | 50 mm, 40 mm, 20 mm    |  |
| Height     | 12.3 cm / 4.84 in | Interpupillary adjustment               | 48 mm - 74 mm          |  |
| Weight     | 239 g / 0.52 lbs  | Ophthalmoscope filters                  | Clear, red-free, amber |  |

To request a quote or a demo, speak with your local sales representative or email us at retina@iridex.com.





**Emergo Europe**Prinsessegracht 20
2514 AP The Hague
The Netherlands



IRIDEX • 1212 Terra Bella Avenue • Mountain View, CA 94043 • U.S.A. 800.388.4747 (U.S.A. only) • 1.650.940.4700 • retina@iridex.com • iridex.com/lio

© 2020 IRIDEX. All rights reserved. IRIDEX, the IRIDEX logo, OcuLight, IQ 577, IQ 532 and trademarks or registered trademarks of IRIDEX. LT0731.A 12.2020





<sup>\*</sup> May vary with refractive power. Retinal spot when treating an emmetropic patient eye with 20 D lens.

<sup>\*\*</sup> At aerial image plane.

<sup>&</sup>lt;sup>1</sup>The LIO Plus is indicated for transpupillary retinal photocoagulation.

<sup>&</sup>lt;sup>2</sup> Hasanreisoglu M, Saktanasate J, Schwendeman R, Shields JA, Shields CL. Indocyanine green-enhanced transpupillary thermotherapy for retinoblastoma: Analysis of 42 tumors J Pediatr Ophthalmol Strabismus, 2015;52(6):348-54.

<sup>&</sup>lt;sup>3</sup> Shields CL, Santos MC, Diniz W, Gunduz K, Mercado G, Cater JR, Shields JA. Thermotherapy for retinoblastoma Arch Ophthalmol, 1999;117(7):885-93.