

# Wall mounted instrument column

## IC-1 Series



# The IC-1 saves you valuable floor space

## IC-1

This IC-1 instrument column is wall mounted, resulting in an exceptionally small footprint, saving you valuable floor space and making it easy access for cleaning purpose. The IC-1 is available in two colours and can be set up for one or even two instruments. The vertical elevation of the tabletop allows for seated or standing operation of your instruments.

With the additional option of a monitor and keyboard support, this IC-1 instrument column will turn into a complete flexible workstation, which also allows easy wheelchair access. The IC-1 in combination with Topcon's joystick free instruments is the perfect solution to provide complete refractive eye care on a single square meter.

- » Modern design
- » Tabletop available for one & two instruments
- » Space saving
- » Simple & easy operation of the tabletop
- » Easy wheelchair accessibility
- » Vertical elevation of tabletop by 600mm
- » Standing & sitting position usability
- » Workstation with monitor & keyboard support



The IC-1 wall mounted instrument column is designed especially for combining it with Topcon's line of joystick free, touch screen operated instruments. These instruments can be placed in a conventional setup, against a wall, or even in a corner of a room, allowing

for the operator to be located anywhere around the patient. The modern and innovative design of the IC-1 will seamlessly fit in any consulting room. Complete with touch screen PC and keyboard your IC-1 will turn into a complete workstation.



**red dot award 2014**  
winner

The IC-1 received the Red Dot award 2104 which is an internationally well-known design award.

## EASY CLEANING

Being wall mounted, the IC-1 provides no obstruction for cleaning the floor in a professional way.



## MONITOR/KEYBOARD SUPPORT

The monitor can be placed on either side of the column in variable height positions, the keyboard tray can be rotated and be stored underneath the tabletop.



## COLOUR CONFIGURATIONS

The IC-1 is available in various color settings.

### POLE



### SIDE COVERS



### TABLE TOP



## FOOTPRINT

Saving valuable floorspace in your practice or shop with the extremely small footprint of the IC-1. The installation of even two instruments will take no more than 0,5 m<sup>2</sup> (780mm x 650mm)



## MULTIPLE INSTRUMENTS

The IC-1 is available with a choice of tabletops suitable for one or two instruments.



## VARIOUS INSTRUMENT POSITIONS

The IC-1 tabletop can be operated from three sides. This allows you to install any joystick operated instrument in a "classic" way with patient and operator on each of the sides.



## ELEVATED TABLETOP

The tabletop can be elevated by 600mm to allow the operator to test at sitting or standing position with ease. The minimum tabletop height of only 720mm allows easy examination of children. The up/down button can be installed on any tabletop side of your choice



## WHEELCHAIR ACCESS

Free from the typical table footprint, the clear floor space makes wheelchair patients manoeuvring and positioning very convenient.



# IC-1E, great solution for small spaces

## IC-1E

IC-1E, Instrument column, is part of the successful IC-1 series. In comparison to the IC-1, the IC-1E enables you to use the complete tabletop as a work area. The IC-1E is available with a tabletop for one instrument or a large tabletop which accommodates two instruments. The tabletop comes in two colors: white & blue.

The IC-1E is unique in its compact dimension and small footprint. Similar to the IC-1, the IC-1E, Instrument column, has a remarkable large stroke of 600mm. IC-1E is especially designed with Topcon's joystick free touchscreen instruments in mind.

- » Innovative unique design
- » Extremely space saving
- » Easy wheelchair accessibility
- » Large tabletop elevation stroke of 600mm
- » Standing and sitting position usability
- » Tabletop available for one or two instruments



## FEATURES

## IC-1E

### COLOR CONFIGURATIONS

The IC-1E is available with two different tabletop colors.


### POLE



### TABLE TOP




### IC-1E Dimensions overview

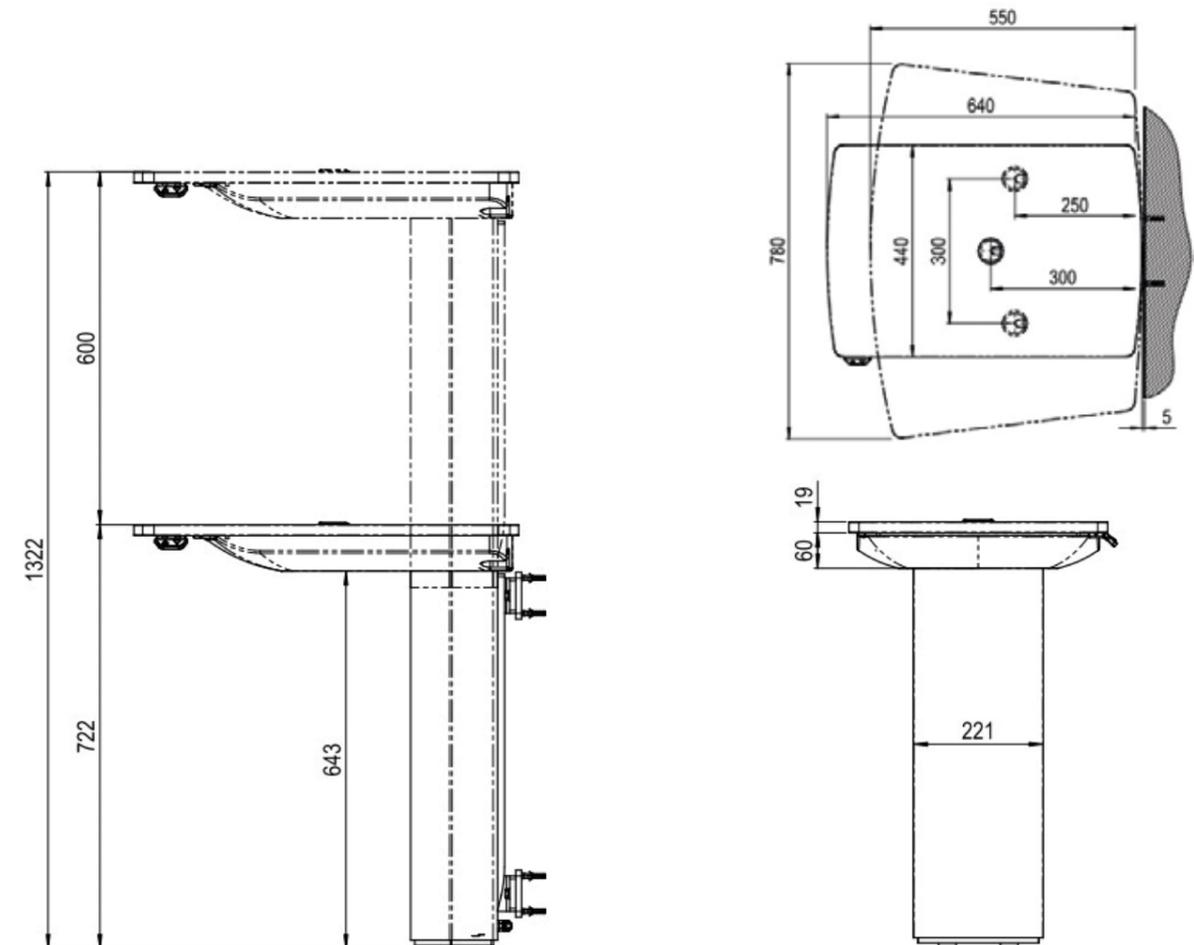

## Specifications

	IC-1	IC-1E
<b>Total power supply</b>	220-240 V / 50-60 Hz	90-240 V / 48-63 Hz
<b>Maximum power consumption</b>	900VA @ 230V	1300VA @ 230V / 650VA @ 110V
<b>Electrical classification</b>	Class I, type B	Class I, Type B
<b>IP class</b>	IP21	IP21
<b>Raising stroke</b>	600 mm	600mm
<b>Maximum load two instrument tabletop</b>	40 kg	40 Kg
<b>Net / Gross weight</b>	53/68 kg	35/48 Kg
<b>Working cycle of the motor</b>	1 / 9 minutes p/h	1 / 6 minutes p/h
<b>Plug connection</b>	Mandatory Earth connection	Mandatory Earth connection

## IC-1 Dimensions overview

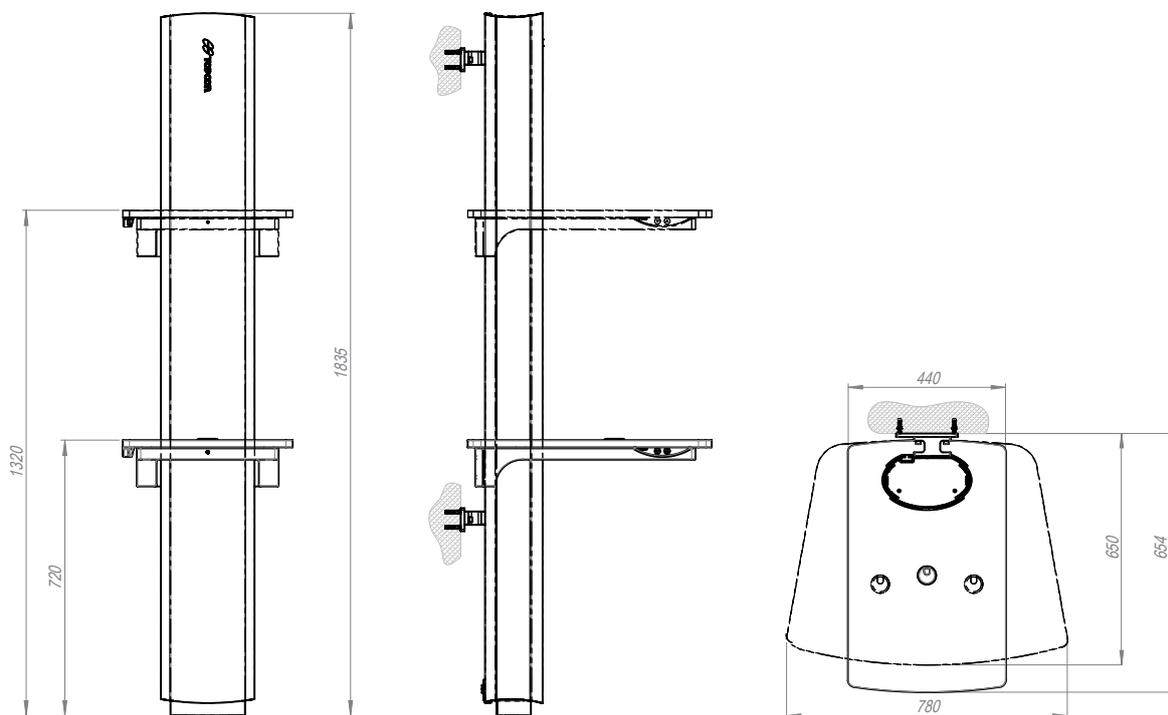

red dot award 2014  
winner

### IMPORTANT

Subject to change in design and/or specifications without advanced notice.

In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation.

#### Topcon Europe Medical B.V.

Essebaan 11; 2908 LJ Capelle a/d IJssel; P.O. Box 145;  
2900 AC Capelle a/d IJssel; The Netherlands  
T: +31-(0)10-4585077; F: +31-(0)10-4585045  
medical@topcon.eu; www.topcon-medical.eu

#### Topcon Danmark

Praestemarksvej 25; 4000 Roskilde, Danmark  
T: +45-46-327500; F: +45-46-327555  
info@topcon.dk; www.topcon.dk

#### Topcon Scandinavia A.B.

Neogatan 2; P.O. Box 25; 43151 Mölndal, Sweden  
T: +46-(0)31-7109200; F: +46-(0)31-7109249  
medical@topcon.se; www.topcon.se

#### Topcon España S.A.

Frederic Mompou, 4;  
08960 Sant Just Desvern; Barcelona, Spain  
T: +34-93-4734057; F: +34-93-4733932  
medica@topcon.es; www.topcon.es

#### Topcon Italy

Viale dell' Industria 60;  
20037 Paderno Dugnano, (MI) Italy  
T: +39-02-9186671; info@topcon.it;  
www.topcon.it

#### Topcon France

BAT A1; 3 route de la révolte, 93206 Saint Denis Cedex  
T: +33-(0)1-49212323; F: +33-(0)1-49212324  
topcon@topcon.fr; www.topcon-medical.fr

#### Topcon Deutschland Medical GmbH

Hanns-Martin-Schleyer Strasse 41;  
D-47877 Willich, Germany  
T: (+49) 2154-885-0; F: (+49) 2154-885-177  
Info@topcon-medical.de; www.topcon-medical.de

#### Topcon Polska Sp. z o.o.

ul. Warszawska 23; 42-470 Siewierz; Poland  
T: +48-(0)32-670-50-45; F: +48-(0)32-671-34-05  
www.topcon-polska.pl

#### Topcon (Great Britain) Ltd.

Topcon House; Kennet Side; Bone Lane; Newbury  
Berkshire RG14 5PX; United Kingdom  
T: +44-(0)1635-551120; F: +44-(0)1635-551170  
medical@topcon.co.uk; www.topcon.co.uk

#### Topcon Ireland

Unit 276, Blanchardstown; Corporate Park 2  
Ballycoolin; Dublin 15, Ireland  
T: +353-18975900; F: +353-18293915  
medical@topcon.ie; www.topcon.ie



### TOPCON EUROPE MEDICAL B.V.

Essebaan 11, 2908 LJ Capelle a/d IJssel, The Netherlands.  
Phone: +31-(0)10-4585077, Fax: +31-(0)10-4585045, E-mail: medical@topcon.eu, www.topcon-medical.eu